

TheReporter

Reprinted from Volume 15, No. 4, 1996

T296034

© 1999 Sigma-Aldrich Co.

For more information, or current prices, contact your nearest Supelco subsidiary listed below. To obtain further contact information, visit our website (www.sigma-aldrich.com), see the Supelco catalog, or contact Supelco, Bellefonte, PA 16823-0048 USA.

ARGENTINA • Sigma-Aldrich de Argentina, S.A. • Av. Pueyrredon 2446/50 • Piso 5-B • Buenos Aires 1119
AUSTRALIA • Sigma-Aldrich Pty. Ltd. • Unit #2, 14 Anella Avenue • Castle Hill NSW 2154
AUSTRIA • Sigma-Aldrich Handels GmbH • Hebbelplatz 7 • A-1110 Wien
BELGIUM • Sigma-Aldrich N.V./S.A. • K. Cardijnplein 8 • B-2880 Bornem
BRAZIL • Sigma-Aldrich Quimica Brasil Ltda. • Rua Sabara, 566-Conj. 53 • 01239-010 São Paulo, SP
CANADA • Sigma-Aldrich Canada, Ltd. • 2149 Winston Park Dr., Oakville, ON L6H 6J8
CZECH REPUBLIC • Sigma-Aldrich s.r.o. • Pobrezni 46 • 186 21 Praha 8
DENMARK • Sigma-Aldrich Denmark A/S • Vejlegaardsvej 65B • DK-2665 Vallensbaek Strand
FINLAND • Sigma-Aldrich Finland/YA-Kemia Oy • Teerisuonkuja 4 • FIN-00700 Helsinki
FRANCE • Sigma-Aldrich Chimie • Chromatographie Supelco • L'Isle d'Abeau Chesnes - B.P. 701 • 38297 Saint-Quentin Fallavier Cedex
GERMANY • Sigma-Aldrich Chemie GmbH • Geschäftsbereich Supelco • Grünwalder Weg 30 • D-82041 Deisenhofen
GREECE • Sigma-Aldrich (o.m.) Ltd. • 72 Argonafton Str. • 16346 Ilioupoli, Athens
HUNGARY • Sigma-Aldrich Kft. • Nagy Diófa u. 7., IV fl. • H-1067 Budapest
INDIA • Sigma-Aldrich Co. • Survey No. 31/1, Sitharamapalaya • Mahadevapura P.O. • Bangalore 560 048
IRELAND • Sigma-Aldrich Ireland Ltd. • Airton Road • Tallaght • Dublin 24
ISRAEL • Sigma Israel Chemicals Ltd. • Park Rabin • Rohovot 76100
ITALY • Sigma-Aldrich s.r.l. • Via Gallarate, 154 • 20151 Milano
JAPAN • Sigma-Aldrich Japan K.K. • Division Supelco • J.L. Nihonbashi Building • 1-10-15 Nihonbashi Horidome-cho, Chuo-ku • Tokyo 103
KOREA • Sigma-Aldrich Korea • Samhan Camus Annex, 10th Floor • 17-26 Yoido-dong, Yungdeungpo-ku • Seoul
MALAYSIA • Sigma-Aldrich (M) Sdn. Bhd. • 9-2, Jalan 2/128, Taman Gembira • Off Jalan Kuchai Lama • 58200 Kuala Lumpur • Selangor
MEXICO • Sigma-Aldrich Quimica S.A. de C.V. • Calle 6 North No. 107 • Parque Industrial Toluca 2000 • 50200 Toluca
NETHERLANDS • Sigma-Aldrich Chemie BV • Postbus 27 • 3330 AA Zwijndrecht
NORWAY • Sigma-Aldrich Norway • Sandakerveien 102 • N-0483 Oslo
POLAND • Sigma-Aldrich Sp. z o.o. • Szelagowska 30 • 61-626 Poznań
PORTUGAL • Sigma-Aldrich Quimica, S.A. • P.O. Box 131 • Sintra 2710
RUSSIA • Sigma-Aldrich Russia • TOO Techmedbiochem • Makarenko Str. 2/21 • Building 1, Flat 22 • Moscow 103062
SINGAPORE • Sigma-Aldrich Pte. Ltd. • 102E Pasir Panjang Road • #08-01 Citilink Warehouse • Singapore 118529
SOUTH AFRICA • Sigma-Aldrich (pty) Ltd. • CNR Kelly & Ackerman Streets • Southern Life Industrial Park Unit • Unit 16/17 • Jet Park 1459
SPAIN • Sigma-Aldrich Quimica, S.A. • Apt. Correos 161 • 28100 Alcobendas, Madrid
SWEDEN • Sigma-Aldrich Sweden AB • Solkraftsvägen 14C • 135 70 Stockholm
SWITZERLAND • Supelco Switzerland • Industriestrasse 25 • P.O. Box 260 • CH-9471 Buchs
UNITED KINGDOM • Sigma-Aldrich Company Ltd. • Supelco UK • Fancy Road, Poole • Dorset BH12 4QH
UNITED STATES • Supelco • Supelco Park • Bellefonte, PA 16823-0048 • Phone 800-247-6628 or 814-359-3441 • Fax 800-447-3044 or 814-359-3044 • email:supelco@sial.com

H

 **SUPELCO**

ISO 9001
REGISTERED

Supelco is a member of the Sigma-Aldrich family. Supelco products are sold through Sigma-Aldrich, Inc. Sigma-Aldrich warrants that its products conform to the information contained in this and other Sigma-Aldrich publications. Purchaser must determine the suitability of the product for a particular use. Additional terms and conditions may apply. Please see the reverse side of the invoice or packing slip.


Analyze Stack Emissions, Using VOST Air Sampling Tubes and Capillary GC

I. DeGraff

United States Environmental Protection Agency SW-846, Method 0030, Volatile Organic Sampling Train (VOST), describes the sampling of emissions from hazardous waste incinerator stacks. Supelco's VOST air sampling tubes are designed to fulfill the method requirements and provide low background levels.

VOST air sampling tubes (Figure A) meet the specifications for US EPA SW-846, Method 0030, Volatile Organic Sampling Train, which explains how to sample for volatile compounds (boiling point range 30-100°C) in emissions from hazardous waste incinerator stacks. The VOST tubes — VOST 100 and VOST 200 — trap these emissions for thermal desorption and gas chromatography/mass spectrometry (GC/MS) analysis.

Figure A. Structure of a VOST Tube

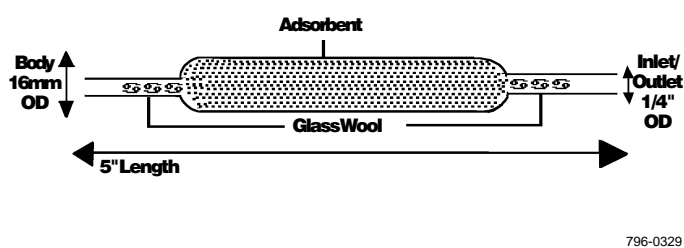

The high-purity adsorbents used in VOST tubes provide consistently low background levels (Table 1). VOST 100 tubes are packed with Tenax® TA 35/60, a porous material based on 2,6-diphenylene oxide polymer that traps volatile and semivolatile compounds and has a low affinity for water and methanol. VOST 200 tubes are packed with Tenax TA 35/60 and petroleum charcoal 20/40 (2:1 by volume). Petroleum charcoal is a naturally occurring adsorbent that traps the more volatile compounds.

The two tubes are assembled in a sampling train in series, with the VOST 100 followed by the VOST 200 as a back-up to trap the more volatile compounds. At flow rates between 250-1000mL/min, a total volume of 5-20L of air is collected on each VOST pair. Up to six VOST pairs can be used to collect samples from the same source, by replacing each pair after a maximum of 20L has been collected.

Table 1. Quality Assurance Specifications

Background:	<20ng single component by GC/FID
Surrogate Recovery:	50-150%
Pressure Drop:	<8 inches water at 50cc/min

Following sample collection, the VOST tubes are thermally desorbed by heating and purging with an inert gas. The effluent is transferred and bubbled into organic-free water in a purge-and-trap system. Analysis is by GC/MS as described in Method 5041, *Protocol for Analysis of Sorbent Cartridges from VOST*. This analysis technique is used because the sampled tubes might contain high quantities of moisture, which could otherwise present problems in a thermal desorption system. More than the 35 compounds listed in Method 5041 can readily be collected with the VOST tubes. However, low recoveries of polar compounds may result because of poor purging efficiency from water.

Figure B compares a calibration standard at 20ng with blank VOST tubes that have been spiked with 50ng chlorobenzene-d₅ as a surrogate to monitor extraction efficiency. The chromatograms illustrate the high purity of the adsorbent, preconditioned and ready to use. The GC conditions encompass the elution range of the target compounds from the EPA method. Peaks comprising the calibration standard represent the most commonly found adsorbent contaminants.


Each VOST tube is sealed with stainless steel Swagelok® end fittings and enclosed in a screw-capped glass storage container that holds a charcoal packet to maintain sample integrity. The storage container is then packed in a protective flex tube. Each batch of VOST tubes is tested for recovery, background, and pressure drop. Quality specifications are listed in Table 1. In addition, each tube has a certificate of analysis and a unique serial number, which ensures accurate sample identification. Tubes can be reused following conditioning with inert gas at temperatures 10-15°C above the desorption temperature.

We can make VOST tubes to your specifications. Choose your own type of tube material, dimensions, or adsorbent. Quotations and prototypes are available upon request. Contact our Technical Service chemists for more information.

Reprinted from Volume 15, No. 4, 1996


Figure B. VOST Tubes Exhibit Low Background

Sampling Tubes: **VOST 100 (Cat. No. 20074-U)**
VOST 200 (Cat. No. 20075-U)
 Purge Trap: **VOCARB™-4000**
 Cat. No.: **20308**
 Column: **SPB™-624, 75m x 0.53mm ID x 3.0µm film**
 Cat. No.: **25432**
 Oven: 35°C (5min) to 200°C at 10°C/min, hold 2min
 Carrier: helium, 8mL/min
 Det.: FID, 300°C

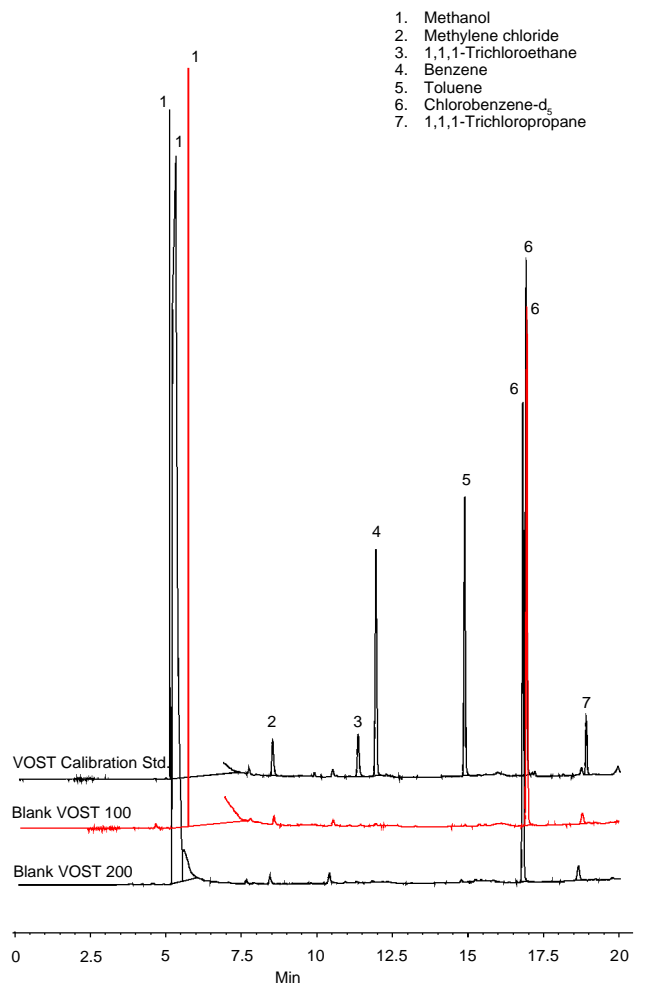

796-0314, 0315, 0316

Ordering Information:

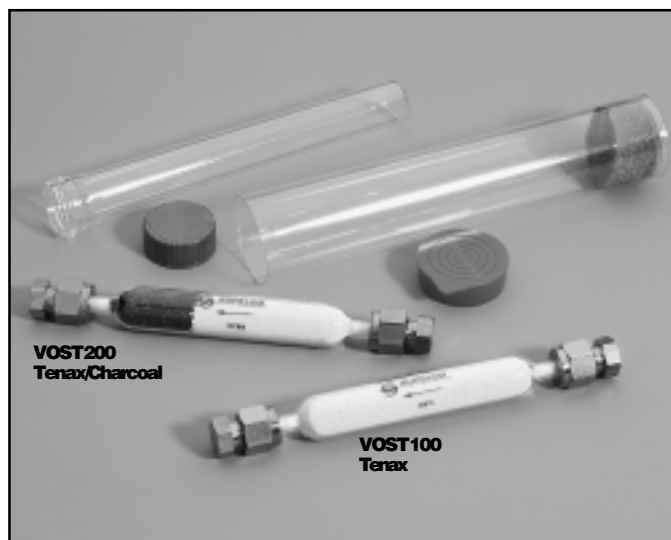

996-0072

Description	Cat. No.
VOST 100 (Tenax TA)	20074-U
VOST 200 (Tenax TA and petroleum charcoal)	20075-U
Empty Glass Tube, 11.5cm x 16mm ID	21993
Glass Storage Container	21998
SPB-624 Fused Silica Capillary Column, 75m x 0.53mm ID, 3.0µm film	25432
VOCARB-4000 Purge Trap	20308*

*For Tekmar® LSC-1, LSC-2, LSC-2000, LSC-4000 instruments. Refer to our catalog for traps for other instruments.

Trademarks

Swagelok — Crawford Fitting Co.
 Tekmar — Tekmar Co.
 Tenax — Enka Research Institute Arnhem
 SPB, VOCARB — Sigma-Aldrich Co.

Fused silica columns manufactured under HP US Pat. No. 4,293,415.



Reprinted from Volume 15, No. 4, 1996