

TheReporter

Reprinted from Volume 16, No. 5, 1997

T297025

© 1999 Sigma-Aldrich Co.

For more information, or current prices, contact your nearest Supelco subsidiary listed below. To obtain further contact information, visit our website (www.sigma-aldrich.com), see the Supelco catalog, or contact Supelco, Bellefonte, PA 16823-0048 USA.

ARGENTINA • Sigma-Aldrich de Argentina, S.A. • Av. Pueyrredon 2446/50 • Piso 5-B • Buenos Aires 1119
AUSTRALIA • Sigma-Aldrich Pty. Ltd. • Unit #2, 14 Anella Avenue • Castle Hill NSW 2154
AUSTRIA • Sigma-Aldrich Handels GmbH • Hebbelplatz 7 • A-1110 Wien
BELGIUM • Sigma-Aldrich N.V./S.A. • K. Cardijnplein 8 • B-2880 Bornem
BRAZIL • Sigma-Aldrich Quimica Brasil Ltda. • Rua Sabara, 566-Conj. 53 • 01239-010 São Paulo, SP
CANADA • Sigma-Aldrich Canada, Ltd. • 2149 Winston Park Dr., Oakville, ON L6H 6J8
CZECH REPUBLIC • Sigma-Aldrich s.r.o. • Pobrezni 46 • 186 21 Praha 8
DENMARK • Sigma-Aldrich Denmark A/S • Vejlegaardsvej 65B • DK-2665 Vallengbaek Strand
FINLAND • Sigma-Aldrich Finland/YA-Kemia Oy • Teerisuonkuja 4 • FIN-00700 Helsinki
FRANCE • Sigma-Aldrich Chimie • Chromatographie Supelco • L'Isle d'Abeau Chesnes - B.P. 701 • 38297 Saint-Quentin Fallavier Cedex
GERMANY • Sigma-Aldrich Chemie GmbH • Geschäftsbereich Supelco • Grünwalder Weg 30 • D-82041 Deisenhofen
GREECE • Sigma-Aldrich (o.m.) Ltd. • 72 Argonafton Str. • 16346 Ilioupoli, Athens
HUNGARY • Sigma-Aldrich Kft. • Nagy Diófa u. 7., IV fl. • H-1067 Budapest
INDIA • Sigma-Aldrich Co. • Survey No. 31/1, Sitharamapalaya • Mahadevapura P.O. • Bangalore 560 048
IRELAND • Sigma-Aldrich Ireland Ltd. • Airton Road • Tallaght • Dublin 24
ISRAEL • Sigma Israel Chemicals Ltd. • Park Rabin • Rohovot 76100
ITALY • Sigma-Aldrich s.r.l. • Via Gallarate, 154 • 20151 Milano
JAPAN • Sigma-Aldrich Japan K.K. • Division Supelco • JL Nihonbashi Building • 1-10-15 Nihonbashi Horidome-cho, Chuo-ku • Tokyo 103
KOREA • Sigma-Aldrich Korea • Samhan Camus Annex, 10th Floor • 17-26 Yoido-dong, Yungdeungpo-ku • Seoul
MALAYSIA • Sigma-Aldrich (M) Sdn. Bhd. • 9-2, Jalan 2/128, Taman Gembira • Off Jalan Kuchai Lama • 58200 Kuala Lumpur • Selangor
MEXICO • Sigma-Aldrich Quimica S.A. de C.V. • Calle 6 North No. 107 • Parque Industrial Toluca 2000 • 50200 Toluca
NETHERLANDS • Sigma-Aldrich Chemie BV • Postbus 27 • 3330 AA Zwijndrecht
NORWAY • Sigma-Aldrich Norway • Sandakerveien 102 • N-0483 Oslo
POLAND • Sigma-Aldrich Sp. z o.o. • Szelagowska 30 • 61-626 Poznań
PORTUGAL • Sigma-Aldrich Quimica, S.A. • P.O. Box 131 • Sintra 2710
RUSSIA • Sigma-Aldrich Russia • TOO Techmedbiochem • Makarenko Str. 2/21 • Building 1, Flat 22 • Moscow 103062
SINGAPORE • Sigma-Aldrich Pte. Ltd. • 102E Pasir Panjang Road • #08-01 Citilink Warehouse • Singapore 118529
SOUTH AFRICA • Sigma-Aldrich (pty) Ltd. • CNR Kelly & Ackerman Streets • Southern Life Industrial Park Unit • Unit 16/17 • Jet Park 1459
SPAIN • Sigma-Aldrich Quimica, S.A. • Apt. Correos 161 • 28100 Alcobendas, Madrid
SWEDEN • Sigma-Aldrich Sweden AB • Solkraftsvägen 14C • 135 70 Stockholm
SWITZERLAND • Supelco Switzerland • Industriestrasse 25 • P.O. Box 260 • CH-9471 Buchs
UNITED KINGDOM • Sigma-Aldrich Company Ltd. • Supelco UK • Fancy Road, Poole • Dorset BH12 4QH
UNITED STATES • Supelco • Supelco Park • Bellefonte, PA 16823-0048 • Phone 800-247-6628 or 814-359-3441 • Fax 800-447-3044 or 814-359-3044 • email:supelco@sial.com

H

This article is archived from a past issue of The Supelco Reporter. Information in the article was appropriate at the time of publication, but product specifications, catalog numbers, and availability may have changed over time.

If you have questions about applying methodology described in this article to a current application, please contact our technical service chemists.



Supelco is a member of the Sigma-Aldrich family. Supelco products are sold through Sigma-Aldrich, Inc. Sigma-Aldrich warrants that its products conform to the information contained in this and other Sigma-Aldrich publications. Purchaser must determine the suitability of the product for a particular use. Additional terms and conditions may apply. Please see the reverse side of the invoice or packing slip.

An Innovative Approach to Coffee Bean Characterization, Using SEC and Principal Component Analysis

J. Dominicki, International Marketing, Supelco, Bellefonte, PA, USA


Coffee beans are one of the largest traded commodities at international levels. Their quality traditionally has been judged by specially trained individuals who taste, smell, and visually examine the beans (sensory analysis). Although today's international trade requirements are stricter, methods for detecting trading fraud and chemical standards for characterizing coffee beans have yet to be developed. Efficient methods for characterizing coffee beans are becoming a necessity.

Investigators at the Universidade Federal do Rio de Janeiro, Instituto de Química characterized coffee beans by using size exclusion chromatography (SEC), with two forms of detection, and principal component analysis (PCA). SEC excludes, or fractions, compounds according to their molecular weight — molecules with larger molecular weights are eluted before those with smaller molecular weights. PCA is a statistical method that analyzes highly correlated variables using nonlinear multiple regression analysis, partial regression analysis, and least square techniques.

Six samples of green coffee beans from various countries of origin were first identified as high or low quality by sensory analysis — beans A, B, and C were of high quality; beans D, E, and F were of low quality. Each of the six samples was separated on a TSKgel G3000SW column, and analyzed with a UV detector and a differential refraction index detector. The TSKgel G3000SW column was calibrated using protein standards, to enable the investigators to split the chromatograms into several regions by molecular weight. The chromatograms in Figure A present a distinct profile showing the low and the high molecular weight regions for a low-quality bean.

The first region, up to 21 minutes, consists of molecules with larger molecular weights, greater than 1.2×10^4 . The second region, 21 to 40 minutes, consists of molecules with low molecular weights, less than 1.2×10^4 . This region consists of very important substances — total chlorogenic acid,

Figure A. Low-Quality Coffee Bean by SEC

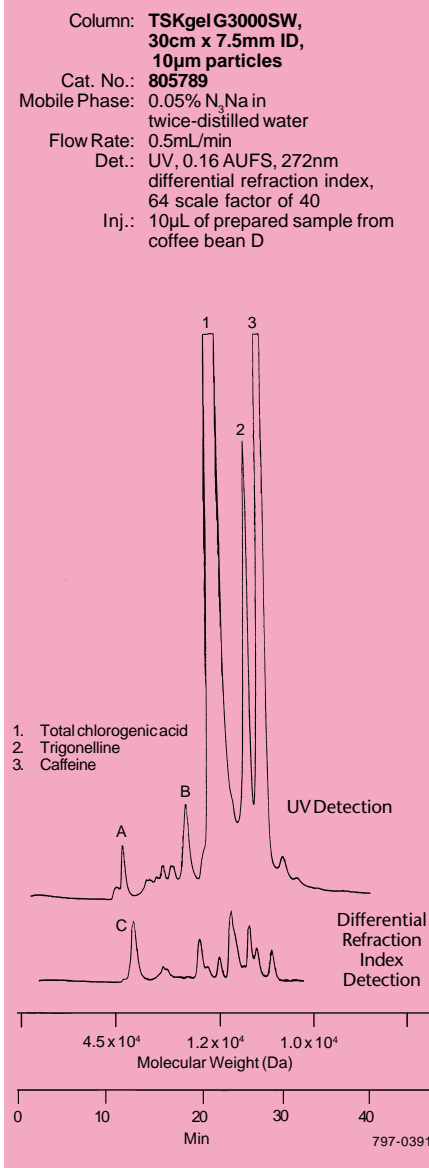

trigonelline, and caffeine — shown respectively in Figure A as peaks 1, 2, and 3 (1).

The compounds shown as peaks A, B, and C contribute 90% of the varieties' quality attributes among the samples — their peak areas are the principal parameters representing the variability of the analyses.

Therefore, these small peaks represent the final quality characteristics of coffee preparations.

Results of the chromatographic analyses were processed by a statistical method — PCA. In PCA three methods are used to assess discrimination stability: principal component analysis, modified PCA, and discriminant function analysis with cross validation. With PCA, variation of a complex system of variables can be expressed by a simple, graphical representation. System integrity and reproducibility is maintained as long as the original sample conditions are kept constant.

Results from PCA showed that SEC, particularly with UV detection, provided a remarkable separation. Furthermore, SEC is capable of simultaneous analysis of the important peaks, whereas reversed-phase HPLC requires different chromatographic conditions for each compound, and headspace-GC analysis results in a high number of crowded peaks.

Further studies are necessary to provide detailed identification of components. Additional complementary analytical techniques, such as mass spectrometry and GC or LC, can help us obtain a detailed characterization profile of most aroma and bouquet compounds for coffee.

Ordering Information:

Description	Cat. No.
TSKgel G3000SW Gel Filtration Columns 30cm x 7.5mm ID, 10 μ m particles	805789

Reference

- De Maria, C. A. B., L. C. Trugo, and R. F. A. Moreira, *Food Chemistry*, 52: 447 (1995). Reference not available from Supelco.

Information in this article used with permission from Sociedade Brasileira de Química, Universidade Federal do Rio de Janeiro, Instituto de Química, Brazil. Professors Ricardo F.A. Moreira, Carlos A.B. De Maria, and Luiz C. Trugo.

TSK-GEL is a registered trademark of Tosoh Corp.