

The Essential Role of Enriched Stable Nuclides in Positron Emission Tomography (PET)



Professor Michael J. Welch
Division of Radiological Sciences
Washington University School of Medicine

Positron Emission Tomography (PET) is one of the most rapidly expanding areas of medical diagnosis. PET is routinely used in determining the extent of tumor metastases, utilizing an analog of glucose in which the hydroxy group at the two position is replaced by a radioactive fluorine atom. This compound, 2-fluoro-2-deoxy-D-glucose (FDG), is taken up in tissues in an amount very similar to glucose, but is retained in the tissue unlike normal glucose, which is rapidly metabolized largely to water and carbon dioxide. This retention of the radioactivity allows the visualization of tumors, which metabolize glucose to a greater extent than normal tissue. In research studies, the glucose utilization of the tumors can be quantified. Due to the two-hour half-life of fluorine-18, this compound must be distributed from regional centers. All of these centers utilize enriched oxygen-18 to produce the fluoride used to make FDG. Fluorine-18 is reacted with agents also supplied by Sigma-Aldrich to produce the radioactive drug used in the medical diagnostic studies.

Other enriched compounds include enriched molecular oxygen and nitrogen enriched in nitrogen-13 and are also used as target materials to produce PET radiopharmaceuticals. Oxygen gas is used in several centers to produce another radioactive drug, fluorine-18-labeled dopamine, used to study brain function, particularly in patients with schizophrenia and other psychiatric diseases. In the most common cyclotron used to produce PET radiopharmaceuticals, the enriched nitrogen is used to produce a series of very simple compounds containing radioactive oxygen-15. Oxygen-15, which has a two-minute half-life, can be used to study brain blood flow and brain oxygen metabolism. In a multicenter study funded by the National Institutes of Neurological Diseases and Stroke, labeled oxygen is being used to predict which patients will benefit from a brain surgery intracranial/extracranial bypass surgery, which restores blood flow to areas of the brain where delivery of oxygen is significantly reduced.

Imaging relies upon the supply of enriched stable isotopes. In the future, it is anticipated that a whole battery of new compounds will be used to study many other parameters in tumors, brain and heart diseases, and they will rely upon the continued supply of stable isotopes.

ISOTECH™ —Your Ultimate Source for Materials for the Production of ¹⁸FDG

- Water-¹⁸O, 97 atom % ¹⁸O (329878)
- ¹⁵N₂/O₂ gas mix ratio 39:1 (600911)
- Mannose Triflate (M1568)
- Kryptofix® (291110)



PET Kits—Ready to Use

Designed for the MX_{FDG} and FX_{FDG} Synthesizer Module

676055 Full Reagent Kit for FX_{FDG} Module

676047 Full Reagent Kit for MX_{FDG} Module

We offer **Custom Packaging** to serve your needs and save you money—ask us how.

Literature of Interest

1. Welch, M. J. et al. Assessment of myocardial metabolism in diabetic rats using small-animal PET: a feasibility study. *J. Nucl. Med.* **2006**, 47(4), 689.
2. Welch, M. J., et al. Production of nonstandard PET radionuclides and the application of radiopharmaceuticals labeled with these nuclides. *Ernst Schering Res Found Workshop*, **2007**, 62, 159. Review.
3. Karantanis, D. et al. Contribution of F-18 FDG PET-CT in the detection of systemic spread of primary central nervous system lymphoma. *Clin. Nucl. Med.*, **2007**, 32(4), 271.
4. Vansteenkiste, J. F.; Stroobants, S. S. PET scan in lung cancer: current recommendations and innovation. *J. Thorac. Oncol.*, **2006**, 1(1), 71. Review.
5. Kubota, K.; Yukihiro, M.; Ito, K. Perspective for imaging of tumor metabolism. *Nippon Rinsho.* **2007**, 65(2), 352. Review. Japanese.
6. Landau, B. R., et al. 6-Fluoro-6-deoxy-D-glucose as a Tracer of Glucose Transport. *Am. J. Physiol. Endocrinol. Metab.*, **2007**, 3.
7. Cherry, S. R. Fundamentals of positron emission tomography and applications in preclinical drug development. *J. Clin. Pharmacol.*, **2001**, 41(5), 482. Review.
8. Johannsen, B. The usefulness of radiotracers to make the body biochemically transparent. *Amino Acids.*, **2005**, 29(4), 307. Review.



PET Applications

Sigma-Aldrich offers target materials for the production of radionuclides ^{18}F , ^{11}C , and ^{15}O . In addition to the target materials, we also offer all reagents used in the production of ^{18}F -FDG. PET research studies have various applications that include Cardiology, Oncology, Neurology, Pharmacology, and Neuropsychology and Cognitive Neuroscience.

Cardiology



PET is the most common form of diagnostic testing for use in cardiology, for example myocardial viability and coronary heart disease. The radionuclide ^{13}N , produced from Water- ^{18}O (**329886**) is used in myocardial perfusion. PET scans help to differentiate viable myocardium from infarcted tissue in patients that are suspected to have hibernating or stunned myocardium.

Oncology



The onset of PET imaging has generated various interests with investigative studies in oncology using ^{18}F -FDG. ISOTEC supplies the Water- ^{18}O (**329878**), Kryptofix® (**291110**), and the Mannose triflate (**M1568**) necessary to produce the ^{18}F -FDG. Today, PET has distinctive uses in clinical areas such as cancers, lymphoma, melanoma, strokes, and tumors.

Neurology



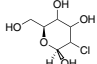
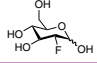
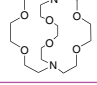
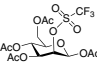
PET is also used to distinguish recurrent brain tumors from radiation fibrosis or necrosis, which is essential to neurological studies. Oxygen- $^{15}\text{O}_2$, produced from $^{15}\text{N}_2$ (99%)/Oxygen (RG) Gas Mix 39:1 (**600911**), measures the flow of blood to the upper or main portion of the brain. The isotope ^{11}C , produced from $^{14}\text{N}_2$ (**608661**), is used as a radiotracer in PET scans to study normal/abnormal brain functions. It is one of the methods in localizing areas of the brain affected by epileptic seizures.

Pharmacology and Neuropsychology



PET is widely used in pre-clinical and clinical trials to study psychiatric and cognitive disorders along with developing new radiolabel drugs.

Compounds of Interest

Structure	Cat. No.	Compound	MW
	C203	2-Chloro-2-deoxy-D-glucose	198.6
	F5006	2-Deoxy-2-fluoro-D-glucose	182.2
	291110	Kryptofix® 4,7,13,16,21,24-Hexaoxa-1,10-diazabicyclo[8.8.8]hexacosane, 98%	376.5
	M1568	Mannose Triflate β-D-Mannopyranose 1,3,4,6-tetra-O-acetate 2-O-trifluoromethanesulfonate 1,3,4,6-Tetra-O-acetyl-2-O-trifluoromethanesulfonyl-β-D-mannopyranose TATM	480.4
$^{15}\text{N}_2/\text{O}_2$	600911	$^{15}\text{N}_2$ (98%)/ O_2 (RG) gas mix ratio 39:1	Minimum 98 atom % ^{15}N
$^{18}\text{O}_2$	602892	Oxygen- $^{18}\text{O}_2$ (Gas)	Minimum 99 atom % ^{18}O
H_2^{18}O	329878	Water- ^{18}O	Minimum 97 atom % ^{18}O

Order All Your
PET Products From
ONE SOURCE.

ISOTEC™

For more information visit us at
sigma-aldrich.com/isotec
or E-mail us at isosales@sial.com
Order and Technical Service:
1-800-448-9760 • 1-937-859-1808

World Headquarters

3050 Spruce St., St. Louis, MO 63103
(314) 771-5765
sigma-aldrich.com

Order/Customer Service (800) 325-3010 • Fax (800) 325-5052

Technical Service (800) 325-5832 • sigma-aldrich.com/techservice

Development/Bulk Manufacturing Inquiries SAFC™ (800) 244-1173

The SIGMA-ALDRICH Group



©2007 Sigma-Aldrich Co. All rights reserved.

SIGMA, SAFC, SAFC™, SIGMA-ALDRICH, ISOTEC, ALDRICH, FLUKA, and SUPELCO are trademarks belonging to Sigma-Aldrich Co. and its affiliate Sigma-Aldrich Biotechnology, L.P. Riedel-de Haën® trademark under license from Riedel-de Haën GmbH. Sigma brand products are sold through Sigma-Aldrich, Inc. Sigma-Aldrich, Inc. warrants that its products conform to the information contained in this and other Sigma-Aldrich publications. Purchaser must determine the suitability of the product(s) for their particular use. Additional terms and conditions may apply. Please see reverse side of the invoice or packing slip. Kryptofix is a registered trademark of Merck Schuchardt OHG.



SIGMA-ALDRICH

Accelerating Customers' Success
through Leadership in Life Science,
High Technology and Service

IOA
02859-504866
0057