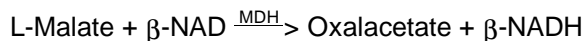


SIGMA QUALITY CONTROL TEST PROCEDURE

**Enzymatic Assay of CITRATE SYNTHASE
(EC 4.1.3.7)****PRINCIPLE:**

Abbreviations used:

MDH = Malic Dehydrogenase

β -NAD = β -Nicotinamide Adenine Dinucleotide, Oxidized Form

β -NADH = β -Nicotinamide Adenine Dinucleotide, Reduced Form

CoA = Coenzyme A

CS = Citrate Synthase

CONDITIONS: T = 37°C, pH = 8.0, A_{340nm}, Light path = 1 cm

METHOD: Continuous Spectrophotometric Rate Determination

REAGENTS:

- A. 100 mM Tris Buffer, pH 8.0 at 37°C
(Prepare 100 ml in deionized water using Trizma Base, Prod. No. T-1503. Adjust to pH 8.0 at 37°C with 1 M HCl.)
- B. 100 mM L-Malic Acid
(Prepare 10 ml in Reagent A using L(-)Malic Acid, Monosodium Salt, Prod. No. M-1125.)
- C. 50 mM β -Nicotinamide Adenine Dinucleotide, Oxidized Form Solution (β -NAD)
(Prepare 2 ml in deionized water using β -Nicotinamide Adenine Dinucleotide, from Yeast, Prod. No. N-7004.)
- D. 2.0 mM Acetyl Coenzyme A Solution (Acetyl CoA)
(Prepare 2.0 ml in Reagent A using Acetyl Coenzyme A, Sodium Salt, Prod. No. A-2056.)
- E. Malic Dehydrogenase Enzyme Solution (MDH)
(Prepare 1.0 ml in Reagent A containing approximately 560 units/ml Malic Dehydrogenase, Prod. No. M-2634.)

**Enzymatic Assay of CITRATE SYNTHASE
(EC 4.1.3.7)**

REAGENTS: (continued)

- F. Citrate Synthase Enzyme Solution (CS)
(Prepare a solution containing 0.2 - 0.4 units of Citrate Synthase in cold Reagent A.)

PROCEDURE:

Prepare a reaction cocktail by pipetting (in milliliters) the following reagents into a suitable container:

Reagent A (Buffer)	7.50
Reagent B (L-Malic Acid)	2.00
Reagent C (β -NAD)	1.00
Reagent D (Acetyl CoA)	2.00
Reagent E (MDH)	1.00
Deionized Water	15.50

Mix and adjust to pH 8.0 at 37°C with 1 M HCl or 1 M NaOH, if necessary.

Pipette (in milliliters) the following reagents into suitable cuvettes:

	<u>Test</u>	<u>Blank</u>
Reaction Cocktail	2.90	2.90

Mix by inversion and equilibrate to 37°C. Monitor the $A_{340\text{nm}}$ until constant using a suitably thermostatted spectrophotometer. Then add:

Reagent F (CS)	0.10	-----
Reagent A (Buffer)	-----	0.10

Immediately mix by inversion and record the increase in $A_{340\text{nm}}$ for approximately 5 minutes. Obtain the $\Delta A_{340\text{nm}}/\text{minute}$ using the maximum linear rate for both the Test and Blank.

CALCULATIONS:

$$\text{Units/ml enzyme} = \frac{(\Delta A_{340\text{nm}}/\text{min Test} - \Delta A_{340\text{nm}}/\text{min Blank})(3)(\text{df})}{(6.22)(0.1)}$$

3 = Volume (in milliliters) of the assay

df = Dilution factor

6.22 = Millimolar extinction coefficient of β -NADH at 340 nm

0.1 = Volume (in milliliter) of enzyme used

**Enzymatic Assay of CITRATE SYNTHASE
(EC 4.1.3.7)**

CALCULATIONS: (continued)

$$\text{Units/mg solid} = \frac{\text{units/ml enzyme}}{\text{mg solid/ml enzyme}}$$

$$\text{Units/mg protein} = \frac{\text{units/ml enzyme}}{\text{mg protein/ml enzyme}}$$

UNIT DEFINITION:

One unit will form 1.0 μmole of citrate from oxalacetate and acetyl CoA per minute at pH 8.0 at 37°C.

FINAL ASSAY CONCENTRATION:

In a 3.0 ml reaction mix, the final concentrations are 46 mM Tris, 6.7 mM L-malic acid, 1.7 mM β -nicotinamide adenine dinucleotide, 0.13 mM acetyl CoA, 56 units malic dehydrogenase, and 0.02 - 0.04 unit citrate synthase.

REFERENCE:

Srere, P.A. (1969), *Methods in Enzymology*, XIII, 1-11

NOTES:

1. This assay is based on the cited reference.
2. Where Sigma Product or Stock numbers are specified, equivalent reagents may be substituted.

Sigma warrants that the above procedure information is currently utilized at Sigma and that all Sigma-Aldrich, Inc. products conform to the information in this and other Sigma-Aldrich, Inc. publications. Purchaser must determine the suitability of the information and product(s) for their particular use. Additional terms and conditions may apply. Please see reverse side of the invoice or packing slip.