

## Product Information

# β-Nicotinamide adenine dinucleotide 2'-phosphate reduced tetrasodium salt hydrate

≥97% (HPLC)

**N7505**

## Product Description

CAS Registry Number: 2646-71-1


Synonyms: β-NADPH, Coenzyme II reduced tetrasodium salt, 2'-NADPH hydrate, NADPH, TPNH, Triphosphopyridine nucleotide reduced tetrasodium salt, NADPH Na<sub>4</sub>, TPNH<sub>2</sub> Na<sub>4</sub>, Dihyronicotinamide adenine dinucleotide phosphate tetrasodium salt

Molecular Formula: C<sub>21</sub>H<sub>26</sub>N<sub>7</sub>Na<sub>4</sub>O<sub>17</sub>P<sub>3</sub> • xH<sub>2</sub>O

Formula Weight: 833.35 (anhydrous basis)

E<sup>m</sup><sub>M</sub> (340nm) = 6.22 (pH >10)A<sub>260nm</sub> / A<sub>340nm</sub> = 2.32

Structure:

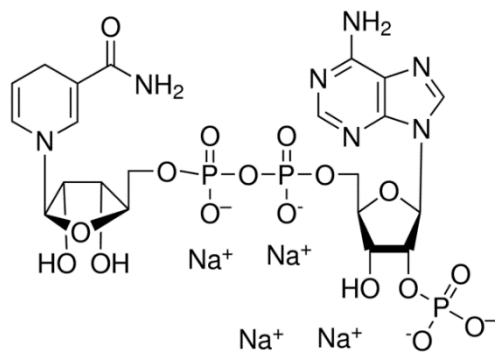

β-Nicotinamide adenine dinucleotide 2'-phosphate (β-NADPH) is a product of the pentose phosphate pathway, a multifunctional pathway whose primary purpose is to generate reducing power, in the form of β-NADPH. β-NADPH transfers H<sup>+</sup> and 2e<sup>-</sup> to oxidized precursors in the reduction reactions of biosynthesis. Thus, β-NADPH cycles between catabolic and biosynthetic reactions, and serves as the carrier of reducing power in the same way that ATP serves as the carrier of energy.<sup>1</sup>

Enzymes that use β-NADPH as a coenzyme include glutathione reductase, diacetyl reductase, dihydrofolate reductase, glutamic dehydrogenase, *p*-hydroxybenzoate hydroxylase, NADPH-FMN oxidoreductase, nitrate reductase and thioredoxin reductase. β-NADPH is also involved with cytochrome P450 electron transport systems.<sup>2</sup>

This listing of β-Nicotinamide adenine dinucleotide phosphate reduced form, N7505, is prepared by the chemical reduction of β-Nicotinamide adenine dinucleotide phosphate. Several theses<sup>3-5</sup> and dissertations<sup>6-14</sup> have cited use of N7505 in their research protocols.

## Precautions and Disclaimer

For R&D use only. Not for drug, household, or other uses. Please consult the Safety Data Sheet for information regarding hazards and safe handling practices.

## Storage/Stability

It is recommended to store N7505 desiccated at -20 °C, protected from light. The normal impurities and/or decomposition products are β-NADP and Monophosphoadenosine 5'-diphosphoribose.

It is suggested to prepare β-NADPH solutions fresh and use promptly, unless you are sure that this is an unnecessary precaution for your work. However, it has been reported that a 0.5 mM solution in 0.02 M NaOH (pH 12.3) showed no loss of purity in a week at 4 °C or -85 °C, but a 13% loss at -20 °C.<sup>15</sup> One publication has investigated the solution stability of NADPH.<sup>16</sup>

## Preparation Instructions

β-NADPH is tested for solubility in 0.01 M sodium hydroxide at 50 mg/mL.

## References

1. Wood, W.B. *et al.*, *Biochemistry: A Problems Approach*. W.A. Benjamin, Inc. (Menlo Park, CA), p. 195 (1974).
2. Berg, J.M. *et al.*, *Biochemistry*, 5<sup>th</sup> edition. W.H. Freeman and Co. (New York, NY), p. 564 (2002).
3. Olsgard, Mandy Lee, "Toxicological evaluation of inhalation exposure to benzene and toluene in a raptorial bird, the American kestrel, *Falco sparverius*". University of Saskatchewan, M.Sc. thesis, pp. 112, 113 (2007).
4. Huang, Huang, "Biological Activities of Selected Dietary *Apiaceae*, *Lamiaceae*, and *Fabaceae*". Carleton University, M.Sc. thesis, p. 16 (2010).
5. Scott, Kirsten Lee, "The Potential Toxic Effects of Chronic Doxorubicin Treatment on the Rat Pancreas and the Role of Ghrelin in this Context". Stellenbosch University, M.Sc. thesis, p. 53 (2018).
6. Mohapatra, Sridev, "Understanding polyamine metabolism through transgenic manipulation in poplar suspension cultures". University of New Hampshire, Ph.D. dissertation, p. 96 (2008).
7. Chrysostomou, Constantine, "Addressing Intrinsic Challenges for Next Generation Sequencing of Immunoglobulin Repertoires". University of Texas at Austin, Ph.D. dissertation, p. 125 (2014).
8. Beyerle, Jolantha, "Assessment of mRNA, protein levels and activities of xenobiotic metabolizing enzymes in colon and rectal mucosa of colorectal cancer patients". Ruperto-Carola University (Heidelberg, Germany), Dr. rer. nat. dissertation, pp. 118, 119 (2016).
9. Gimpfl, Martina Christiane, "The role of fatty acid quality in a maternal obesogenic diet for adipose tissue development in offspring mice". Technischen Universität München, Dr. rer. nat. dissertation, p. 107 (2016).
10. Hord, Jeffrey Mack, "Triggering NOX2 / ROS signaling during disuse-induced skeletal muscle atrophy". Texas A&M University, Ph.D. dissertation, p. 92 (2016).
11. Kuehne, Jennifer, "Anti-inflammatory properties of New Zealand honey". University of Auckland, Ph.D. dissertation, pp. 27, 47 (2016).
12. Litten, Brett Andrew, "Biphasic droplet microfluidics in relation to pharmaceutical industrial biochemical screening". University of Manchester, Ph.D. dissertation, p. 81 (2016).
13. Vassallo, Joanne, "The bioavailability and biological effects of nanomaterials towards *Escherichia coli*: with reference to the soil environment". University of Plymouth, Ph.D. dissertation, p. 46 (2018).
14. Sharick, Joseph Thomas, "Development of Optical Metabolic Imaging for Pancreatic and Breast Cancer Precision Medicine". University of New Hampshire, Ph.D. dissertation, p. 49 (2019).
15. Passonneau, J.V., and Lowry, O.H., *Enzymatic Analysis. A Practical Guide*. Humana Press (Totowa, NJ), p. 15 (1993).
16. Wu, J.T. *et al.*, *Clin. Chem.*, **32(2)**, 314-319 (1986).

## Notice

We provide information and advice to our customers on application technologies and regulatory matters to the best of our knowledge and ability, but without obligation or liability. Existing laws and regulations are to be observed in all cases by our customers. This also applies in respect to any rights of third parties. Our information and advice do not relieve our customers of their own responsibility for checking the suitability of our products for the envisaged purpose.

The information in this document is subject to change without notice and should not be construed as a commitment by the manufacturing or selling entity, or an affiliate. We assume no responsibility for any errors that may appear in this document.

## Technical Assistance

Visit the tech service page at [SigmaAldrich.com/techservice](https://SigmaAldrich.com/techservice).

## Standard Warranty

The applicable warranty for the products listed in this publication may be found at [SigmaAldrich.com/terms](https://SigmaAldrich.com/terms).

## Contact Information

For the location of the office nearest you, go to [SigmaAldrich.com/offices](https://SigmaAldrich.com/offices).

The life science business of Merck KGaA, Darmstadt, Germany operates as MilliporeSigma in the U.S. and Canada.

MilliporeSigma, and Sigma-Aldrich are trademarks of Merck KGaA, Darmstadt, Germany or its affiliates. All other trademarks are the property of their respective owners. Detailed information on trademarks is available via publicly accessible resources.

© 2022 Merck KGaA, Darmstadt, Germany and/or its affiliates. All Rights Reserved.

N7505pis Rev 03/22 DAG,NA,GCY,MAM

**MILLIPORE  
SIGMA**