


Monoclonal Anti- β -Tubulin, Isotype III, Clone SDL-3D10 (T 8660)

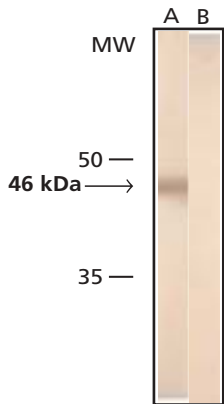

Immunoblot

Whole cell extract of rat brain cytosol was separated on SDS-PAGE and blotted with Monoclonal Anti- β -Tubulin Isotype III (Product No. T 8660). The antibody was developed with Goat anti-Mouse IgG, Peroxidase conjugate (Product No. A 9917) and AEC substrate.

Lane A: Antibody dilution 1:400

Lane B: Negative control (only secondary antibody)