




SIGMA QUALITY CONTROL TEST PROCEDURE

Product Information

Enzymatic Assay of PHOSPHATASE, ACID, PROSTATIC (EC 3.1.3.2)

PRINCIPLE:

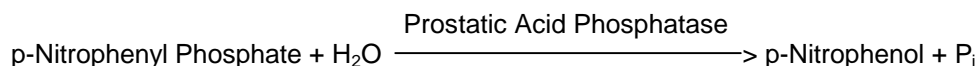

Abbreviation used:

P_i = Inorganic phosphate

CONDITIONS: T = 37°C, pH = 4.8, A_{410nm}, Light path = 1 cm

METHOD: Spectrophotometric Stop Rate Determination

REAGENTS:

- A. 90 mM Citrate Buffer Solution, pH 4.8 at 37°C
(Use Citrate Buffer Solution, Sigma Stock No. 104-4.)
- B. 40 mM Tartrate Buffer, with 90 mM Citrate, pH 4.8 at 37°C.
(Use Tartrate Acid Buffer Solution, Sigma Stock No. 104-12.)¹
- C. 11.3 mM p-Nitrophenyl Phosphate (PNPP)
(Prepare 5 ml in deionized water using p-Nitrophenyl Phosphate Sigma 104 Phosphatase Substrate, Sigma Stock No. 104-0.)
- D. 100 mM Sodium Hydroxide Solution (NaOH)
(Prepare 50 ml in deionized water using Sodium Hydroxide, Anhydrous Sigma Prod. No. S-5881.)
- E. Prostatic Acid Phosphatase Enzyme Solution
(Immediately before use, prepare a solution containing 0.015 - 0.02 unit/ml of Acid Phosphatase in cold deionized water.)

**Enzymatic Assay of PHOSPHATASE, ACID, PROSTATIC
(EC 3.1.3.2)**

PROCEDURE:

Pipette (in milliliters) the following reagents into suitable containers:

	<u>Test1</u>	<u>Test2</u>	<u>Blank</u>
Reagent A (Citrate Buffer)	0.50	-----	0.50
Reagent B (Tartrate Buffer)	-----	0.50	-----
Reagent C (PNPP) 0.50	0.50	0.50	

Mix and equilibrate to 37°C. Then add:

	<u>Test1</u>	<u>Test2</u>	<u>Blank</u>
Reagent E (Enzyme Solution)	0.20	0.20	-----

Immediately mix and incubate at 37°C for exactly 30 minutes. Then add:

Reagent D (NaOH) 3.80	3.80	3.80	
Reagent E (Enzyme Solution)	-----	-----	0.20

Mix and record the A_{410nm} for both the Tests and Blank in a suitable spectrophotometer.

CALCULATIONS:

Total Acid Phosphatase Activity:

$$\text{Units/ml Total enzyme} = \frac{(\Delta A_{410nm}/\text{min Test}_1 - \Delta A_{410nm}/\text{min Blank})(5.0)(df)}{(30)(18.3)(0.2)}$$

Nonprostatic Acid Phosphatase Activity:

Units/ml Nonprostatic enzyme =

$$\frac{(\Delta A_{410nm}/\text{min Test}_2 - \Delta A_{410nm}/\text{min Blank})(5.0)(df)}{(30)(18.3)(0.2)}$$

Prostatic Acid Phosphatase Activity:

Units/ml Prostatic Enzyme = Units/ml Total Enzyme - Units/ml Nonprostatic Enzyme

Enzymatic Assay of PHOSPHATASE, ACID, PROSTATIC (E C 3.1.3.2)

CALCULATIONS: (continued)

5.0 = Total volume (in milliliters) of solution

30 = Time of assay (in minutes)

18.3 = Millimolar extinction coefficient of p-Nitrophenol at 410 nm at an Alkaline pH

0.2 = Volume (in milliliter) of enzyme used

df = Dilution factor

$$\text{Units/mg solid} = \frac{\text{units/ml Prostatic Enzyme}}{\text{mg solid/ml Total Enzyme}}$$

$$\text{Units/mg protein} = \frac{\text{units/ml Prostatic Enzyme}}{\text{mg protein/ml Total Enzyme}}$$

UNIT DEFINITION:

One unit will hydrolyze 1.0 μ mole of p-nitrophenyl phosphate per minute at pH 4.8 at 37°C.

Prostatic acid phosphatase activity is the difference between the total acid phosphatase activity and the acid phosphatase activity in the presence of 20 mM tartrate.

FINAL ASSAY CONCENTRATION:

In a 1.20 ml reaction mix, the final concentrations are 38 mM citric acid, 4.7 mM p-nitrophenyl phosphate and 0.004 unit prostatic acid phosphatase.

NOTES:

1. Tartrate Buffer negates roughly 95% of the prostatic acid phosphatase.
2. Where Sigma Product or Stock numbers are specified, equivalent reagents may be substituted.

Sigma warrants that the above procedure information is currently utilized at Sigma and that all Sigma-Aldrich, Inc. products conform to the information in this and other Sigma-Aldrich, Inc. publications. Purchaser must determine the suitability of the information and product(s) for their particular use. Additional terms and conditions may apply. Please see reverse side of the invoice or packaging slip.